



MULTI-FUNCTION DISPERSION MACHINE LMDM-A10

Multi-function Dispersion Machine LMDM-A10


Multi-function Dispersion Machine LMDM-A10 is 1.2 L stuff jar multi-function, combine function of disperse, mill and mix in one machine. The handing ability ranges from 200 to 8000 mL and offers diameter 300 mm Lift height. Features with digital display and brainpower control. The machine is used in laboratory experiments and for dispersing, grinding and mixing tiny volume materials. The device is equipped with grinding disk and dispersing disk of Step less speed regulation. The container is double-avoidance structure, can be cooled or heated and easy to operate, stable performance, easy to clean and convenient to replace the parts.

FEATURES

- Designed with Multi-function, combine function of disperse, mill and mix in one machine
- Speed ranges from 100 to 8000 rpm
- The disperser head it easy to disassemble and clean
- Equipped with widely batch tackle scope
- Offers digital display and brainpower control

APPLICATIONS

Multi-function Dispersion Machine is used for dispersing, grinding and mixing tiny volume materials across used in Paint, coating, pharmaceutical, food, water treatment, scientific experiments and other industries.

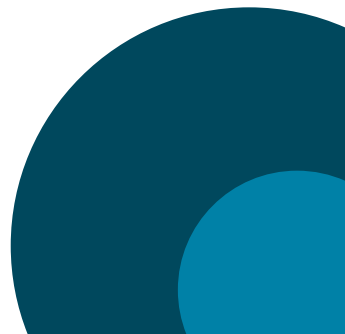

SPECIFICATION

Model	LMDM-A10
Stuff jar	1.2 L
Handing ability	200 to 8000 mL
Speed	100 to 8000 rpm
Material	SS 304
Disc	Disperse: Φ50, Φ60
	Mill: Φ45 Φ60
Lift height	300 mm
Dimension	450 × 500 × 800 mm
Power consumption	400 W
Power supply	220V, 50Hz
Net Weight	25 kg

